



EYEWEAR

Your eyes are 4 times more sensitive to radiation than other parts of your body. Not too long ago, few doctors used leaded eyewear. With the increase in occurrences of cataracts, lead glasses are now required in most procedures using x-ray.

Protech Medical manufactures a wide range of leaded eyewear including: lead glasses, prescription eyewear, goggles, fitovers, and faceshields.

Our dedication to manufacturing in the U.S. and commitment to sourcing the highest quality materials has made us the number one manufacturer in the world.

Eyewear Sizing and Add-Ons:

- Refer to page 61 for eyewear size information and extras.



All lead glasses are manufactured by Protech Medical Inc. in Florida.



Made with industry standard 0.75mm lead equivalent lenses.



Included in all eyewear purchases is a hard case, head strap, and a 1 year manufacturer warranty against defects.



EYEWEAR LINES

- **PREMIUM** - Our premium eyewear lines contain some of the most popular frames we offer. Ranging from our 99 series to our 89, 90, and 42 fitovers. This eyewear line features proven styles and quality frames that have been sold for years.
- **ECONOMY** - We understand leaded eyewear can be costly. So we created a line dedicated to low cost frames combined with Protech Medical's high quality standards.
- **STANDARD** - The standard line consists of all our classic frames. From our standard 53 wrap to the 553 and 662 metals. The standard line is packed with frames that are retro by nature but provide essential lateral protection.
- **BRAND NAME** - We doubt this one needs much introduction but our brand lines are brand name only frames ranging from Nike, Bebe, Gucci, Rayban and Wiley X. These frames have been carefully selected to offer maximum lateral protection while still looking great.

With the accessibility of modern designed frames and ultra-light lenses, protective leaded eyewear is used by roughly 95% of physicians and staff in the following vocations:



CARDIAC CATH



EP LABS



UROLOGY



**INTERVENTIONAL
RADIOLOGY**



PAIN MANAGEMENT



**ORTHOPEDIC
SURGERY**

The typical interventional cardiologist gets exposed to 20 to 30 mSv each year. This level of exposure is nearly 10 times greater than the level of exposure received by the entire body.

- 98 Superlite -

Wraparounds

Premium



Available Colors: Brown

Model: 98
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Material: Durable Nylon
Weight: 63g
Frame Size: 63-15-120



The 98 Superlite is a new lead glasses model. It is designed to offer comfort and provide protection for the entire eye area without the need of an added side shield for splash radiation exposure. The curved front and temples with rubber tips assures a close, secure fit.

- 9935 Ultralite -

Wraparounds

Premium



Rx

Available Colors: Black, Blue Trans., Red Trans., Silver

Model: 9935 / 9935L
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Optional 0.50mm (Regular only)
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range (reg.): +5.00 to -5.00
Prescription Range (lg.): +5.00 to -5.00
Material: Durable Nylon
Weight (reg.): 60g
Weight (lg.): 65g
Frame Size (reg.): 54-21-120
Frame Size (lg.): 58-21-120



BLK



BLU



RD



SIL



The 9935 model is a funky wrap-around frame designed to offer comfort and provide protection for the entire eye area without the need of an added side shield for lateral x-ray exposure. This frame is ideal for those who have smaller facial structures. It features a soft rubber nose bridge providing extra support and increased comfort over long periods of use. The curved front and temples with rubber tips assure a close, secure fit.

- 9941 Ultralite -

Wraparounds


Premium



Rx

Model: 9941 / 994153
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range (reg.): +3.00 to -5.00
Prescription Range (sm.): +6.00 to -9.00
Material: Durable Nylon
Weight (reg.): 64g
Weight (sm.): 60g
Frame Size (reg.): 55-22-130
Frame Size (sm.): 53-22-130

Available Colors: Black, Blue, Red, Silver

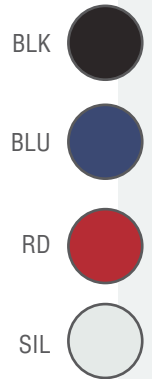

The 9941 Ultralites are our premium wrap-around frame. Designed to offer comfort and provide protection for the entire eye area without the need of added side shield protection, this frame is ideal for those who have a high power prescription. The 9941 also features a soft rubber nose bridge providing extra support for increased comfort and a light-weight feeling.

- 9941 Ultralite w/ Nose Pads -

Wraparounds


Premium



Rx

Model: 9941A
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.00 to -5.00
Material: Durable Nylon
Weight: 64g
Frame Size: 55-22-130

Available Colors: Black, Blue, Red, Silver

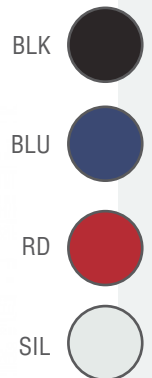

The 9941 Ultralites are our premium wrap-around frame. Designed to offer comfort and provide protection for the entire eye area without the need of added side shield protection, this frame is ideal for those who have a high power prescription. The 9941 also features a soft rubber nose bridge providing extra support for increased comfort and a light-weight feeling.

- 99 Alumilite -

Wraparounds

Premium

CE

Rx

Model: 99AL
 CE Marked: Yes
 Special: Flexible Spring Hinge Temple
 Front Lens: 0.75mm Lead Glass
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +2.50 to -3.00
 Material: Aluminum
 Weight: 70g
 Frame Size: 62-17-150

Available Colors: Black, Cobalt Blue, Gold, Red, Silver



Our model 99 alumilite is constructed out of lightweight aluminum & combined with adjustable non-slip nose pads to provide maximum comfort for all users. The glasses curved shape, flexible spring hinge temples, along with rubber tip temples assure a tight fit and provide splash radiation protection.

- 99 Ultralite -

Wraparounds


Premium

CE

Rx

Model: 99 / 99SM
 CE Marked: Yes
 Front Lens: 0.75mm Lead Glass
 Optional 0.50mm (small only)
 Side Shields: Optional 0.50mm Pb Strips
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range (reg.): +3.00 to -3.00
 Prescription Range (sm.): +3.00 to -5.00
 Material: Durable Nylon
 Weight (reg.): 67g
 Weight (sm.): 59.5g
 Frame Size (reg.): 64-18-140
 Frame Size (sm.): 56-18-140

Available Colors: Black, Blue, Taupe/Brown, Silver, Red

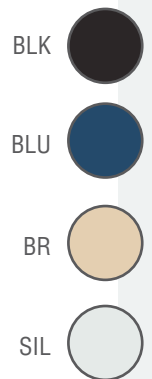

The 99 Ultralites are our premium wrap-around model. They are designed to offer comfort and provide radiation protection for the entire eye area without the need of added side shields for lateral radiation exposure. The 99 has an optional soft rubber nose bridge which can provide extra support and increased comfort. The 99 lead glasses are a staple eyewear at Protech and are one of the most light-weight large sized models.

- 99 Ultralite w/ Nose Bridge -

Wraparounds


Premium



Rx

Available Colors: Black, Blue, Taupe/Brown, Silver, Red

Model: 99R
 CE Marked: Yes
 Front Lens: 0.75mm Lead Glass
 Optional 0.50mm (small only)
 Side Shields: Optional 0.50mm Pb Strips
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range (reg.): +3.00 to -3.00
 Prescription Range (sm.): +3.00 to -5.00
 Material: Durable Nylon
 Weight (reg.): 67g
 Weight (sm.): 59.5g
 Frame Size (reg.): 64-18-140
 Frame Size (sm.): 56-18-140

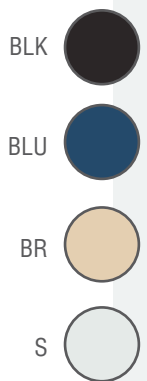

The 99 Ultralites are our premium wrap-around model. They are designed to offer comfort and provide radiation protection for the entire eye area without the need of added side shields for lateral radiation exposure. The 99 has an optional soft rubber nose bridge which can provide extra support and increased comfort. The 99 lead glasses are a staple eyewear at Protech and are one of the most light-weight large sized models.

- Women's 99 Ultralite -

Wraparounds


Premium



Rx

Available Colors: Teal, Purple, Light Pink, Light Blue, Fuschia

Model: L99 / L99SM
 CE Marked: Yes
 Front Lens: 0.75mm Lead Glass
 Optional 0.50mm (small only)
 Side Shields: Optional 0.50mm Pb Strips
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range (reg.): +3.00 to -3.00
 Prescription Range (sm.): +3.00 to -5.00
 Material: Durable Nylon
 Weight (reg.): 67g
 Weight (sm.): 59.5g
 Frame Size (reg.): 64-18-140
 Frame Size (sm.): 56-18-140

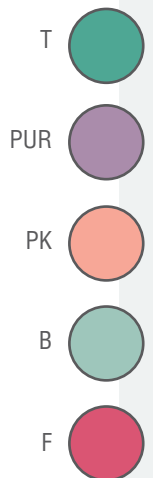

The L99 Ladies Ultralites are our premium wrap-around model for women. They are designed to offer comfort and provide radiation protection for the entire eye area without the need of added side shields. The L99 also features a soft rubber nose bridge, or nose pads which provide extra support and increased comfort. The L99 womens lead glasses are a staple eyewear at Protech and are one of most light-weight radiation eyewear.

- Lattice -

Wraparounds

Designer

Available Colors: Black, Brown

Model: L3273
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Material: Acetate
Weight: 72g
Frame Size: 64-18-135



Women's Lattice lead glasses are part of our new designer line. They feature a classy buckle design on stylish patterned temples. Though the lens shape is rectangular, we have classified these eyewear as a wrap-around frame due to sufficient lateral x-ray protection.

- Plaid -

Wraparounds

Designer

Available Colors: Silver, Brown

Model: L3280
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Material: Acetate
Weight: 68g
Frame Size: 60-17-135



The Women's Plaid frame is so named for its stylish plaid accents. Though the lens shape is rectangular, we have classified these eyewear as a wrap-around frame due to sufficient lateral x-ray protection.


- Zebra -

Wraparounds

Designer

Available Colors: Black, Brown

Model: L4656
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Material: Acetate
Weight: 71g
Frame Size: 61-15-120

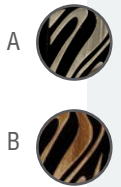

With modern elegance the zebra coloring and patterned temples offer an excellent flavor as designer leaded eyewear. Our Women's Zebra frame is part of our new designer line. Though the lens shape is rectangular, we have classified these eyewear as a wrap-around frame due to sufficient lateral x-ray protection.

- Vader -

Wraparounds

Economy

Rx

Available Colors: Gray, Matte Black, Shiny Black, Clear Brown

Model: NY117
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal
Prescription Range: +3.00 to -3.00
Material: Durable Plastic
Weight: 64g
Frame Size: 64-17-125



The Vaders model is the most flexible wrap-around frame available on the radiation safety market. The Vader's ultra-flex technology provides premium comfort at a fraction of the cost of traditional wrap-around leaded eyewear. As with all wrap-around models, the Vaders provide protection for the entire eye area without the need of added side shields for lateral exposure. Additionally, the Vaders are one of the lightest radiation safety frames given their large area of protection.

- Modern Metal -

Metals

Economy

Rx

Model: RE-002
CE Marked: No
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +6.00 to -9.00
Material: Metal
Weight: 73g
Frame Size: 48-20-135

Available Colors: Gold, Light Gun Metal



The model 002 is our wrap styled metal frame. Its modern design has made it the most popular metal frame. The model 002 features standard lead strip side shields but is also available without them. This frame is the flagship of our economy metal line.

- Metal Aviator -

Metals

Economy

Rx

Model: RE-003-005
CE Marked: No
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +6.00 to -9.00
Material: Metal
Weight: 82g
Frame Size: 54-17-150

Available Colors: Gold



The model 003 is a bold edition to our economy led eyewear lineup. These frames have style and come in our modern gold color. They are also part of the economy line which means they are stylish and affordable!

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- Petite Metal -

Metals

Economy

Rx

Model: RE-007
CE Marked: No
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +6.00 to -9.00
Material: Metal
Weight: 63g
Frame Size: 52-22-150

Available Colors: Gold, Light Gun Metal



005



011



The model 007 is our petite metal economy frame. It features side shields with lead strips and comes in two available colors.

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- Classic Metal -

Metals

Economy

Rx

Model: RE-008
CE Marked: No
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +6.00 to -9.00
Material: Metal
Weight: 64g
Frame Size: 54-19-140

Available Colors: Gold, Nano Blue



005



014



The model 008 is part of our economy line. It comes standard with custom fitted 0.50mm leaded side shields that create a seal to the front lenses for lateral x-ray protection. Being apart of our economy line, this classic metal frame is one of our most affordable lead glasses option. Also available without side shields.

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- 553 Metalite -

Metals

Standard

CE

Rx

Available Colors: Antique Bronze

Model:	553 553S
CE Marked:	Yes
Front Lens:	0.75mm Lead Glass Optional 0.50mm Pb Lens
Side Shields:	Optional 0.50mm Pb Strips
Available Prescription:	Single Vision, Bifocal, Progressive
Prescription Range:	+6.00 to -9.00
Material:	Optical Steel
Weight:	56g
Weight (w/ Side Shields):	62g
Frame Size:	50-18-140



The model 553 is a metal frame featuring custom fitted 0.50mm lead strip side shields that create a seal to the front lenses for necessary protection from oblique angle exposure. Its bronze color and retro style offer a classic look.

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- 557 Metalite -

Metals

Standard

CE

Rx

Available Colors: Gold

Model:	557S
CE Marked:	Yes
Front Lens:	0.75mm Lead Glass Optional 0.50mm Pb Lens
Side Shields:	Optional 0.50mm Pb Strips
Available Prescription:	Single Vision, Bifocal, Progressive
Prescription Range:	+6.00 to -9.00
Material:	Optical Steel
Weight:	63g
Frame Size:	48-20-135



Our model 557S is an optical steel frame featuring custom fitted 0.50mm lead side shields that create a seal to the front lenses for necessary protection from oblique angle exposure to x-ray radiation. The gold color gives this retro style an updated look.

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- 661 Metalite -

Metals

Standard

CE

Rx

Available Colors: Gold

Model: 553
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +6.00 to -9.00
Material: Optical Steel
Weight: 54g
Frame Size: 52-18-144



Our model 661 is a modern styled optical steel frame that comes with optional side shields (662S) as seen below.

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- 662 Metalite -

Metals

Standard

CE

Rx

Available Colors: Gold

Model: 662S
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Optional 0.50mm Pb Lens
Side Shields: Optional 0.50mm Pb Strips
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +6.00 to -9.00
Material: Optical Steel
Weight: 70g
Frame Size: 52-18-145



The model 662S is our most popular metal frame. It features custom fitted side shields and a stylish rectangular frame.

*Includes adjustable non-slip silicon nose pads for maximum comfort.

- 42 Fitover -

Fitovers

Premium



Available Colors: Black

Model: 42BLK
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Side Shields: 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Material: Durable Nylon
Weight: 104g
Frame Size: 63-14-135



The new 42 lead fitover, like the model 89 and 90, is made to wear over prescription eyewear. The 42 Fitover features a smaller vertical dimension than other fitover styles as well as 0.5mm Pb lead sheeting for side shields, thus decreasing the weight by approx. 20%. The fitovers are an excellent choice for those requiring splash protection.

- 89 Medium Fitover -

Fitovers

Premium



Available Colors: Black, Tortoise

Model: 89
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Side Shields: 0.75mm Lead Glass
Available Prescription: None
Prescription Range: None
Material: Durable Nylon
Weight: 107g
Frame Size: 63-11-125



BLK



TORT



Originally designed to fit over prescription eyewear, these leaded fitovers are an excellent choice for those requiring good splash protection from radiation while wearing their normal prescription glasses. The large, clear side shields protect the user from lateral exposure while eliminating the "tunnel" effects of look-alikes. The model 89 is a fraction smaller than the 90 and, consequently, weighs less.

- 90 Large Fitover -


Fitovers

Premium



Available Colors: Black

Model: 90
 CE Marked: Yes
 Front Lens: 0.75mm Lead Glass
 Side Shields: 0.75mm Lead Glass
 Available Prescription: None
 Prescription Range: None
 Material: Durable Nylon
 Weight: 112g
 Frame Size: 66-14-120

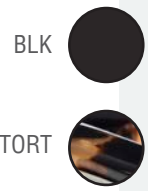

Originally designed to fit over prescription eyewear, model 90 leaded fitovers are an excellent choice for those requiring good splash protection from radiation. The large, clear side shields protect the user from lateral exposure while eliminating the “tunnel” effects of look-alikes. The 90 black is by far our most popular fitover model.

- 53 Wrap -

Classic


Standard



Rx

Available Colors: Tortoise, Red Tortoise, Green Tortoise, Blue Tortoise, Black Tortoise

Model: 53
 CE Marked: Yes
 Special: Nose Pads 53A (Only in tortoise)
 Front Lens: 0.75mm Lead Glass
 Side Shields: 0.75mm Lead Glass
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +6.00 to -9.00
 Material: Durable Nylon
 Weight: 84g
 Frame Size: 58-16-135

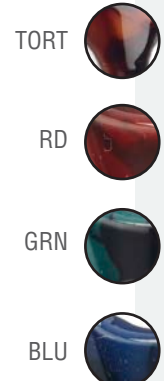

The 53 wrap is constructed of extremely durable nylon that exhibits superior flexibility while minimizing potential breakage. The molded wrap provides maximum comfort with a full saddle bridge to distribute weight evenly over the nose. The lead glass side shields ensure great lateral protection while also providing some peripheral vision.

*Also available with adjustable nose pads (Tortoise only) (53A)

- Bebe Free Spirit -

Designer Brands

Bebe

Rx

Available Colors: Black, Crystal

Model: BB-SPIRIT-7087
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Durable Plastic
Weight: 84g
Frame Size: 58-15-135



001



971



A distinct twist to a classic look, this unique frame comes in two beautiful colors.

- Bebe Fairest -

Designer Brands

Bebe

Rx

Available Colors: Black, Topaz
Crystal, Navy Crystal

Model: BB-FAIR-7086
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Durable Nylon
Weight: 74g
Frame Size: 57-17-130



001



210



455



The Fairest frame offers a unique pairing of color and designwork.

- Bebe Fiery -


Designer Brands

Bebe

Rx

Available Colors: Black, Topaz

Model: BB-FIERY-7082
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Durable Nylon
Weight: 67g
Frame Size: 56-15-125

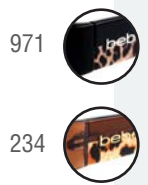

Modern and sexy, the Fiery frame blends a sleek design with luxurious leopard print accents.

- Gucci Wrap -

Designer Brands

Gucci

Rx

Available Colors: Black, Tortoise/
Black

Model: RE-GG-1000S
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Durable Plastic
Weight: 72.5g
Frame Size: 57-15-140



807



WRR



The Gucci Wrap Lead Glasses are a unisex frame featuring patterned accents, a large field of view, and adequate lateral protection from scatter beam radiation.

- Gucci Women's -

Designer Brands

Gucci

Available Colors: Black, Tortoise

Model: RE-GG-3166
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: None
Prescription Range: None
Material: Durable Nylon
Weight: 86g
Frame Size: 57-15-140



D28



791



The Gucci model 3166 is our women's Gucci frame. It features the world renowned Gucci logo on the temple bar with the patented signature web detail.

- Nike Brazen -

Designer Brands

Nike

CE

Rx

Model: N-BR
 CE Marked: Yes
 Front Lens: 0.75mm Lead Glass
 Side Shields: Optional 0.50mm Pb Strips
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -3.50
 Material: Durable Nylon
 Weight: 75g
 Frame Size: 60-16-140

Available Colors: Black w/ Silver,
 Tortoise, Matte Black
 w/ Silver, Mercury Grey w/ Lime,
 Matte Black w/ Blue



571-001

571-204

572-095

758-003

758-049



The Brazen features an advanced ventilated nose bridge for grip, comfort, and reduced fogging. Its lightweight nylon design combines streamlined style with enhanced comfort. The premium Nike Brazen frame offers versatility with its superior coverage and protection.

- Nike Gaze 2 -

Designer Brands

Nike

Model: N-GAZE
 CE Marked: No
 Front Lens: 0.75mm LE
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -3.50
 Material: Lightweight Nylon
 Weight:
 Frame Size:

Available Colors: Matte Crystal
 Dark Grey, Crystal Violet, Black



760-007

836-250

836-001



The Gaze 2 is both sporty and stylish with its cat-eye shape and performance features. Soft rubber inlay nose pads provide enhanced comfort and security, while the frame's six-base design offers a relaxed sporty feel. Secure Wrap Temple technology features special grip zones for better overall fit, comfort, and stability.

- Nike Mercurial 6.0 -

Designer Brands


Nike

CE

Rx

Model: N-MER6
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Lightweight Nylon
Weight: 71g
Frame Size: 61-13-130

Available Colors: Matte Black/
Silver, Midnight Navy, Matte
Crystal Anthracite

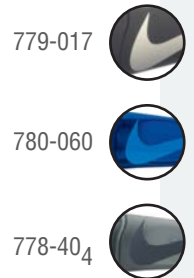

The Mercurial 6.0 is a sporty wrap frame with a bold Nike swoosh on the temples. It features soft rubber inlay nose pads and temple tips for enhanced comfort and grip. The Mercurial's lightweight frame design boasts both sport wrap and secure wrap technology: Maximum coverage and increased peripheral vision combined with grip zones for a better overall fit, comfort, and stability.


- Nike Mercurial 8.0 -

Designer Brands

Nike

Model: N-MER8
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: None
Prescription Range: None
Material: Lightweight Nylon
Weight: 71g
Frame Size: 61-14-135

Available Colors: Pewter Orange,
Dark Royal Blue, White/Black

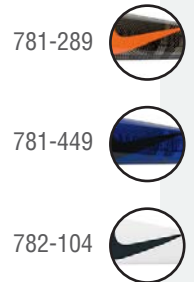

The Mercurial 8.0 is a sporty eight-base frame with a bold Nike swoosh on the temples. It features soft rubber inlay nose pads and temple tips for enhanced comfort and grip. The Mercurial's lightweight frame design boasts both sport wrap and secure wrap technology: Maximum coverage and increased peripheral vision combined with grip zones for a better overall fit, comfort, and stability.

- Nike Premier 6.0 -

Designer Brands

Nike

Rx

Model: N-PRMR
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Lightweight Nylon
Weight: 75g
Frame Size: 59-13-135

Available Colors: Crystal Mercury
Gray, Matte Black/Purple



790-016



791-056



The Premier is a six-based sport frame providing a relaxed, comfortable fit. It features soft rubber inlay nose pads and temple tips for enhanced comfort, grip, and stability. The Premier's Sport Wrap Technology provides maximum coverage and increased peripheral vision.

- Nike Premier 8.0 -

Designer Brands

Nike

Model: N-PRMR8
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: None
Prescription Range: None
Material: Lightweight Nylon
Weight: 75g
Frame Size: 61-14-145

Available Colors: Black/Silver,
Matte Black/Red



792-009



794-065



The Premier is an eight-based sport frame providing a relaxed, comfortable fit with refined curvature. It features soft rubber inlay nose pads and temple tips for enhanced comfort, grip, and stability. The Premier's Sport Wrap Technology provides maximum coverage and increased peripheral vision.

- Nike Poise -

Designer Brands




Nike

Rx

Available Colors: Crystal Gray, Black, Tortoise, Crystal Violet

Model: N-POISE
 CE Marked: No
 Front Lens: 0.75mm Lead Glass
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -3.50
 Material: Lightweight Nylon
 Weight: 65g
 Frame Size: 57-16-135



- 885-001 
- 741-001 
- 741-202 
- 741-205 

The Poise is a modern style with high fashion appeal. Its soft rubber inlay nose pads improve grip and eliminate slippage while its lightweight nylon design provides comfort. The Poise frame features Secure Wrap Technology with special grip zones for a better overall fit, comfort, and stability.

- Nike Rabid -




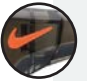
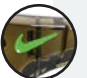
Designer Brands

Nike

Available Colors: Anthracite, Mattle Black/White, Matte Black/Purple, Pewter/Orange

Model: N-RABID
 CE Marked: No
 Front Lens: 0.75mm Lead Glass
 Side Shields: None
 Available Prescription: None
 Prescription Range: None
 Material: Lightweight Nylon
 Weight: 63g
 Frame Size: 63-14-145



- 603-065 
- 604-095 
- 795-056 
- 603-095 
- 603-303 

The Rabid is a high-curvature sports frame providing superior wrap-around coverage and protection. Its light-weight nylon frame provides clarity at all angles during any procedure. Secure Wrap Technology features grip zones throughout the frame for a better overall fit, comfort, and stability.

- Nike Siren -

Designer Brands

Nike

CE

Rx

Model: N-SI
 CE Marked: Yes
 Front Lens: 0.75mm Lead Glass
 Side Shields: Optional 0.50mm Pb Strips
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -3.50
 Material: Durable Nylon
 Weight: 63g
 Frame Size: 58-14-130

Available Colors: Black w/ Black, Tortoise, Deep Royal/Turquoise, Grand Purple/Pink



The Siren frame is designed with style in mind and comes with many features. The Siren's light curves give off a modern look while a ventilated nose bridge provides a high level of comfort and reduces fogging. The Nike Siren's main feature is Secure Wrap Technology with special grip zones providing a better overall fit, comfort, and stability.

- Nike Tour -

Designer Brands

Nike

Model: N-TOUR-744
 CE Marked: No
 Front Lens: 0.75mm Lead Glass
 Side Shields: None
 Available Prescription: None
 Prescription Range: None
 Material: Metal
 Weight: 50g
 Frame Size: 61-13-150

Available Colors: Walnut, Black, Chrome/Blue



The Nike Tour is a high quality sports frame with precision curves. Made of optical quality steel the Tour frame features an adjustable nose bridge for increased comfort. With excellent coverage and enhanced peripheral vision the Tour is premium, light-weight option.

- Nike Varsity -


Designer Brands

Nike

Rx

Model: N-VARS-821
 CE Marked: No
 Front Lens: 0.75mm Lead Glass
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -3.50
 Material: Lightweight Nylon
 Weight: 50g
 Frame Size: 55-14-125

Available Colors: Crystal Blue/
 Yellow Strike, Black/Volt, Matte
 Crystal Gray

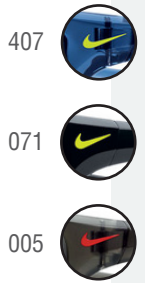

The Varsity is a petite style with sleek curves. Its light-weight, sporty design provides superior coverage and protection.

- Nike Unrest -


Designer Brands

Nike

Rx

Model: N-UNREST
 CE Marked: No
 Front Lens: 0.75mm LE
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -3.50
 Material: Lightweight Nylon
 Weight:
 Frame Size:

Available Colors: Matte Blue
 Pewter, Dark Obsidian White,
 Grey/Red, Black Deep Pewter

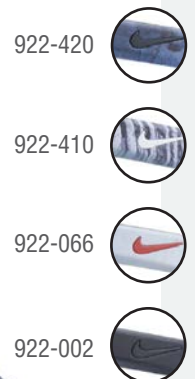

Paying homage to Nike's hometown of Portland, OR, the Unrest frame features unique graphics on the outside temples that are taken from Burnside Skatepark. The graphics are "Special Edition" and digitally printed on the temples.

- Nike 7085 -

Designer Brands

Nike

Rx

Available Colors: Black/Red, Satin Black/Grey, Black/Grey

Model: N-7085
 CE Marked: No
 Front Lens: 0.75mm LE
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -5.50
 Material: Lightweight Nylon
 Weight:
 Frame Size:



24739-001



24739-005



29826-014



The 7085 is a light-weight sports frame built for comfort and durability. It features a ventilated nose bridge for improved comfort and reduced fogging. An adjustable rubber temple sleeve with Nike Lunarlon pattern temple grips provides enhanced stability and added comfort.

- Nike 7087 -

Designer Brands

Nike

Rx

Available Colors: Cargo Khaki/Orange, Obsidian/Game Royal

Model: N-7087
 CE Marked: No
 Front Lens: 0.75mm LE
 Side Shields: None
 Available Prescription: Single Vision, Bifocal, Progressive
 Prescription Range: +3.50 to -5.50
 Material: Lightweight Nylon
 Weight:
 Frame Size:



29826-083



29826-421



The 7087 is a stylish, light-weight frame built for comfort and durability. It features a ventilated nose bridge for improved comfort and reduced fogging. An adjustable rubber temple sleeve with Nike Lunarlon pattern temple grips provides enhanced stability and added comfort.

- Nike 7108 -

Designer Brands

Nike

Rx

Available Colors: Satin Black/Grey,
Satin Crystal Night, Satin
Crystal Teal

Model: N-7108
CE Marked: No
Front Lens: 0.75mm LE
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +5.50 to -7.50
Material: Lightweight Nylon
Weight:
Frame Size:



025



038



320



The 7108 is a premium frame made of ultra light-weight carbon fiber and nylon. While providing excellent comfort and flexibility the 7108 is also durable. Adjustable rubber temple sleeves provide stability and added grip. The 7108 is an optical frame designed for high prescriptions and visual acuity.

- Wiley-X P17 -


Designer Brands

Wiley-X

Rx

Available Colors: Black, Brown

Model: 69696
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: Single Vision, Bifocal, Progressive
Prescription Range: +3.50 to -3.50
Material: Durable Nylon
Weight: 68g
Frame Size: 60-18-125

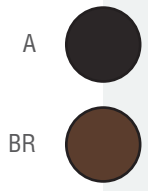

The Wiley X P-17 is similar to the Zak but comes in two colors and is a smaller frame. It also offers wrap around lateral protection from radiation thus eliminating the need for side shields.

- Wiley-X Zak -

Designer Brands

Wiley-X

Available Colors: Black

Model: 65552
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: None
Prescription Range: None
Material: Durable Nylon
Weight: 74g
Frame Size: 65-15-125



The Wiley X Zak is a wide wrap around frame for maximum vision and protection from radiation. Its modern look has contributed to its great popularity and success.

- 68 Goggle -

Goggles

Economy

Available Colors: White/Black

Model: 68
CE Marked: No
Front Lens: 0.75mm Lead Glass
Side Shields: None
Available Prescription: None
Prescription Range: None
Material: Vinyl Polycarbonate
Weight: 102g
Frame Size: 55-19-Adj.



The model 68 is made of soft vinyl that provides a contoured facial fit. They also come with optional lead strip side shields for lateral protection. The mask can fit over petite eyewear in some cases.

- 69 Goggle -

Goggles

Economy



Available Colors: White/Black

Model: 69
CE Marked: Yes
Front Lens: 0.75mm Lead Glass
Side Shields: 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Material: Vinyl Polycarbonate
Weight: 105g
Frame Size: 55-19-Adj.



The model 69 comes standard with lead strip side shields. The frame is made of a soft vinyl that contours to your face.

- 70 Astroflex -

Goggles

Economy

CE

Rx

Available Colors: Clear with Blue

Model:	70/70S
CE Marked:	Yes
Front Lens:	0.75mm Lead Glass
Side Shields:	Optional 0.50mm Pb Strips
Available Prescription:	Single Vision, Bifocal, Progressive
Prescription Range:	+4.00 to -7.00
Material:	Propionate
Weight:	85g
Weight (w/ Side Shields):	88g
Frame Size:	54-19-Adj.



The model 70 is a panoramic lead goggle constructed of impact resistant hypoallergenic plastic propionate. The frame front has molded side shielding and is fitted to the brow for increased radiation splash protection. The temple length and frame angle are adjustable for a uniform custom fit. Model 70S comes with leaded side shields for necessary protection from oblique angle exposure.

- Full Face Shield -

Face Shields



Model:	450FPM
CE Marked:	Yes
Front Shield:	0.10mm Lead Acrylic
Side Shields:	None
Available Prescription:	None
Prescription Range:	None
Material:	Plastic
Weight:	450g
Frame Size:	Adjustable Head Piece



The Full Face Shield uses a lead impregnated protective acrylic shield suspended from a dual adjustable headpiece. The Full face shield provides maximum x-ray protection due to its larger lead impregnated shield. Because of the multi-adjustable headpiece, you can angle the mask to protect from radiation in many different positions as well as adjust the size to comfortably fit your head.

*Mask should be worn with leaded eyewear protection.

- Panoramic Face Shield -

Face Shields



Model:	400PSM
CE Marked:	Yes
Front Shield:	0.10mm Lead Acrylic
Side Shields:	None
Available Prescription:	None
Prescription Range:	None
Material:	Plastic
Weight:	350g
Frame Size:	Adjustable Head Piece



The Panoramic Face Shield utilizes a lead impregnated acrylic shield suspended from a dual adjustable headpiece. This lead face shield provides additional x-ray protection around the sides of the head while sacrificing some front protection. Because of the multi-adjustable headpiece, you can angle the mask to protect from radiation in many different positions as well as adjust the size to comfortably fit your head.

*Mask should be worn with leaded eyewear protection.

EYEWEAR EXTRAS

TEMPLE

The length of your frame's temples is important to ensure stability and comfort.

TEMPLE GRIPS

Available with several eyewear styles, temple grips can provide comfort and stability for the wearer.

LEAD GLASS

Lead glass lenses provide standard 0.75mm Pb Equivalent protection.



NOSE BRIDGE

A well-designed nose bridge is necessary for a proper fit.

NOSE PADS

Ventilated nose pads keep the frame comfortably gripped on your face.

EYEWEAR ACCESSORIES

ITEM	CODE
Case	RE-CASE
Large Case	RE-LCASE
Croakie	RE-CROAK
A-Fog Cleaner	RE-CLEAN
Cloth Wipe	RE-CLOTH

ADD ONS

ITEM	CODE
A/R Coating	RE-AR
Lead Side Shields	RE-SS

FRAME SAMPLES

ITEM	CODE
Non Brand	RE-FS-NB
Nike	RE-FS-N
Bebe	RE-FS-BB
Gucci	RE-FS-G



GLOVES

Protech Medical is the largest manufacturer and distributor of radiation attenuation gloves in the world. We pioneered the development of surgical gloves and lead the way in educating healthcare professionals of their importance. We are proud to offer two types of surgical gloves, the Proguard RR Series and the Proguard 1 & 2. In addition to our surgical gloves we have multiple patient positioning and veterinary style gloves available in various protection options and colors.

A study by the National Center for Biotechnology Information concluded that:

A healthcare worker could receive more than twice the annual limit of radiation exposure simply by holding the head of a trauma patient four times per week with unshielded hands during cervical spine radiography.

The National Council on Radiation Protection & Measurements' annual limit of occupational radiation exposure to: skin, hands, and feet is 500 milli-sievert.

 Rigorously tested

 FDA Approved

- Proguard RR Series-

Attenuation Glove

Exam


CE

FDA

Model: RR1/RR2/RR3
 CE Marked: Yes
 FDA 510k: Yes
 Sterilized: Yes - Ethilene Oxide
 Resterilization: Up to 3 times by Ethilene Oxide
 Allergen Info: Powder Free & Latex Antigen Free
 Features: Anatomically Correct Finger Curvature
 Usage: Procedures requiring high tactile sensitivity

Sizing Info:

- 6 - RR160 RR260
- 6.5 - RR165 RR265
- 7 - RR170 RR270 RR370
- 7.5 - RR175 RR275 RR375
- 8 - RR180 RR280 RR380
- 8.5 - RR185 RR285 RR385
- 9 - RR190 RR290 RR390
- 10 - RR310

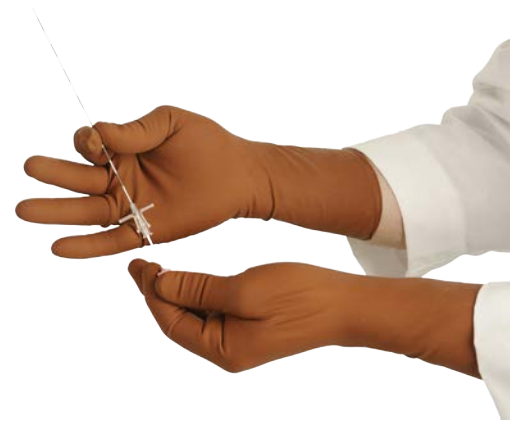

The Proguard radiation attenuation glove is the dominant glove on the market. It is FDA 510k approved, CE certified, and sold worldwide. Used as an exam glove, the Proguard is used in: Radiology, Cath/EP Lab, Orthopedics, Endoscopy, Urological Procedures, Special Procedures, and Pain Management. Attenuates up to 50% more than competitors.

ATTENUATION

	60 kVp	80 kVp	100 kVp	120 kVp	Thickness	Lead Equivalence
Proguard RR1	45%	35%	26%	23%	0.22mm ± 0.02mm	0.016mm Pb
Proguard RR2	55%	43%	35%	31%	0.30mm ± 0.02mm	0.022mm Pb
Proguard RR3	81%	69%	62%	59%	0.60mm ± 0.02mm	0.032mm Pb

BIOCOMPATIBILITY

Test	Standards	Results
Cutaneous Irritation	0 = No erythema/edema 4 = Severe erythema/edema	Index: 0.0
Allergic Sensation	0 = No erythema/edema 4 = Severe erythema/edema	Index: 0.0
Protein Residuals	Modified Lowry	< 50

Due to the presence of natural rubber latex and lead oxide, PROGUARD gloves have undergone tests on cutaneous irritation, allergic sensitization and proteins content with the results shown.

FDA TESTING & EUROPEAN TECHNICAL STANDARDS

Test	Standards
Check for pinholes	ASTM 3578, EN 420
Dexterity	EN 420/03
Protection against ionizing radiations and radioactive contamination	EN 421/94
Protection against chemicals	EN 374/03

- Proguard RR Sleeve -

Attenuation Sleeve

Exam



FDA

Model: SLEEVE
 CE Marked: Yes
 FDA 510k: Yes
 Allergen Info: Powder Free & Latex Antigen Free
 Features: Inner Nitrile coating-reducing proteins, antigens, and irritation
 Sizing Info: Small = SLEEVE/S
 Medium = SLEEVE/M
 Large = SLEEVE/L



The Proguard radiation reducing sleeves are FDA 510(k) approved and may be used in addition to our Proguard gloves. It is meant for ensuring complete x-ray protection from your fingers to your biceps.

ATTENUATION

	60 kVp	80 kVp	100 kVp	120 kVp	Thickness	Lead Equivalence
Proguard Sleeve	53%	43%	35%	31%	0.30mm ± 0.02mm	0.022mm Pb

- Proguard 1 & 2 -

Attenuation Glove

Surgical



FDA



Model: PG-1, PG-2
 CE Marked: Yes
 FDA 510k: Yes
 Sterilized: Yes
 Allergen Info: Powder, Latex, and Lead Free
 Features: Anatomically Correct Finger Curvature
 Usage: Procedures requiring high tactile sensitivity

Sizing Info:

- 6 - PG160 PG260
- 6.5 - PG165 PG265
- 7 - PG170 PG270
- 7.5 - PG175 PG275
- 8 - PG180 PG280
- 8.5 - PG185 PG285
- 9 - PG290 PG290



The Proguard 1 & 2 are lead and latex free.

ATTENUATION

	50 kVp	80 kVp	100 kVp	Thickness	Lead Equivalence
Proguard 1	69%	35%	26%	0.22mm ± 0.02mm	0.22mm Pb
Proguard 2	55%	43%	35%	0.30mm ± 0.02mm	0.22mm Pb

COMING SOON

**BIOCOMPATIBILITY
 &
 FDA TESTING & EUROPEAN TECHNICAL STANDARDS**

- Angio Mitt -

Positioning Mitt



Model:	AP-ANGIOMIT
CE Marked:	Yes
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	0.50, 0.35, or 0.25 LE
Material Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery or Logo
Protection:	0.50mm Lead Equivalence (LE)
Features:	Can be used for left or right hand
Sizing Info:	One Size fits all



Our Angio Mitt is part of a wide selection of radiation shielding accessories. The angio mitt is a flexible mitt used to shield the hand, thus reducing exposure under direct x-ray. This mit can be used on either the right or left hand and is mainly used in positioning.

- Slit Mitt -

Positioning Mitt



Model:	AP-SLITMIT(L/R)
CE Marked:	Yes
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	0.50, 0.35, or 0.25 LE
Material Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery or Logo
Protection:	0.50mm Lead Equivalence (LE)
Features:	Must specify as left or right hand
Sizing Info:	One Size fits all



Our slit mitt is part of a wide selection of radiation shielding accessories. This flexible glove has a lead protective sheet for the top part of your hand and a slit on the bottom to allow for gripping.

- Angiographic Lead Gloves -

Positioning Mitt

Model:	ANGIO
CE Marked:	No
Material:	Leatherette
Customizations:	Embroidery or Logo
Protection:	0.50mm Lead Equivalence (LE)
Features:	Can be used for left or right hand
Sizing Info:	One Size fits all



The Angiographic Lead Glove is mainly used in patient positioning but is also ideal for high radiation exposure procedures. The glove is soft and flexible and can easily be wiped clean.

BARRIERS



Many of the mobile barriers on the market were dated. They also greatly restricted viewing due to large opaque panels. Protech Medical pioneered a new method of producing barriers. First, instead of glass, we used light-weight acrylic. We then made sure that the barrier was composed mostly of acrylic as opposed to a panel. This new, clean looking design is light weight and provides an excellent field of vision.

In line with our commitment to quality, all Protech barriers are produced in the USA. We are now growing our line to also include glass for those in need of additional radiation protection.

- The base features .50mm Pb lead sheeting



0.50mm Lightweight Lead Acrylic



Excellent viewing field

- Full Acrylic Barrier -

Barriers

Mobile

Model:	MB FULL
CE Marked:	No
Protection (Acrylic):	Standard 0.50mm Lead Equivalence
Protection (Base):	Standard 0.50mm Lead Equivalence
Window Dimensions:	55" H x 28" W
Full Dimensions:	72" H x 28" W



Diagnostic x-ray procedures frequently require working personnel to remain in the examination room where they are subject to exposure from secondary radiation. The clear leaded acrylic barrier can provide total-body x-ray shielding during extended procedures. Our clear acrylic barrier rests on an opaque lower panel with smooth-rolling casters giving the barrier great mobility. Its fully clear window gives you a modern alternative to old fashioned shielding barriers. The acrylic lead barrier is shipped unassembled but can be easily assembled using a screwdriver and wrench. For heavy duty radiation protection barriers please see our lead glass barriers.

- Half Acrylic Barrier -

Barriers

Mobile

Model:	MB HALF
CE Marked:	No
Protection (Acrylic):	Standard 0.50mm Lead Equivalence
Protection (Base):	Standard 0.50mm Lead Equivalence
Window Dimensions:	24" H x 28" W
Full Dimensions:	72" H x 28" W



Diagnostic x-ray procedures frequently require working personnel to remain in the examination room where they are subject to exposure from secondary radiation. The clear leaded acrylic barrier can provide total-body x-ray shielding during extended procedures. Our half acrylic barrier has a larger opaque panel which reduces overall weight and cost. With smooth-rolling casters giving the barrier great mobility and a fully clear window allowing excellent viewing, the acrylic barrier is a modern choice for radiation shielding. The acrylic lead barrier is shipped unassembled but can be easily assembled using a screwdriver and wrench. For heavy duty radiation protection barriers please see our lead glass barriers.

- Adjustable Barrier -

Barriers

Mobile

Model:	MB-ADJ
CE Marked:	No
Protection:	Standard 0.50mm Lead Equivalence
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	1.00, 0.75, 0.50, 0.35, or 0.25 LE
Material Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery or Logo
Special:	Custom Sizing Available
Dimensions:	24" W x 48" H
Height Adjustment:	24" to 48" H



Our adjustable mobile radiation barrier can give you maximum shielding from x-ray exposure. Ideal for patient & personnel radiation protection our lead shield is easy to maneuver, and much lighter than the full & half acrylic barriers. This enables greater mobility while maintaining a solid construction with high quality, durable wheel casters.

- Porta Shield -

Barriers

Mobile

Model:	MB-PORTA
CE Marked:	No
Protection:	Standard 0.50mm Lead Equivalence
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	1.00, 0.75, 0.50, 0.35, or 0.25 LE
Material Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery or Logo
Special:	Custom Sizing Available
Dimensions:	24" W x 24" H to 36" H
Height Adjustment:	34" to 60"



The graduated porta shield is a mobile-adjustable lead shield that provides a variety of different adjustable heights blended with a light weight design. The pentagon shaped wheels provide stability and portability while maintaining a clean design. The shield is attached by velcro which enables easy removal for cleaning.

- Graduated Porta Shield -

Barriers

Mobile

Model:	MB-PORTA-GRAD
CE Marked:	No
Protection:	Standard 0.50mm Lead Equivalence
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	1.00, 0.75, 0.50, 0.35, or 0.25 LE
Material Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery or Logo
Special:	Custom Sizing Available
Dimensions:	24" W x 17" H to 36" H
Height Adjustment:	34" to 60" H



The graduated porta shield is similar to the porta shield but differentiates itself with its individual lead flaps. The flaps provide excellent versatility and usage. They are attached by velcro which enable easy removal for cleaning.

- Surgical Barrier -

Barriers

Mobile

Model:	MB-SURGERY
CE Marked:	No
Protection:	Standard 0.50mm Lead Equivalence
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	1.00, 0.75, 0.50, 0.35, or 0.25 LE
Material Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery or Logo
Special:	Custom Sizing Available
Dimensions:	
Height Adjustment:	



The surgical barrier moves with you. It is designed to provide you with maximum protection during procedures requiring fluoroscopy and other sources of x-ray.

- Single/Double-sided Drapes -

Table Accessories

Table Drapes

Available Colors: Refer to pages 27-29.

CE Marked:	Yes
Protection:	0.50mm Lead Equivalence (LE)
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	0.50, 0.35, or 0.25 LE
Fabric Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery & Logo



The single sided table drape uses velcro to adhere to the side of your imaging table. The double sided drape features a saddlebag design for a secure fit to any style imaging table. They provide protection against scatter radiation for your lower body and are specially designed to not obstruct positioning of the C-arm. Each drape attaches to the top by velcro allowing for easy cleaning and adjustability.

	48"x24"	48"x28"	48"x32"	48"x36"	36"x24"	36"x28"	36"x36"
Single sided	TS/48241	TS/48281	TS/48321	TS/48361	TS/36241	TS/36281	TS/3636FP
Double sided	TS/48242	TS/48282	TS/48322	TS/48362	TS/36242	TS/36282	TS/3636DS

- Anesthesia Front Panel -

Table Accessories

Table Drapes

Available Colors: Refer to pages 27-29.

CE Marked:	Yes
Protection:	0.50mm Lead Equivalence (LE)
Core Material Options:	Lead (ST), Light Lead (LW), Lead-Free (LF), Premier (PR)
Core Material Thickness:	0.50, 0.35, or 0.25 LE
Fabric Options:	Nylon, Ripstop, Designer, Vinyl
Customizations:	Embroidery & Logo

Picture Not Available

The anesthesia front panel is not a drape it is a full panel designed to securely fit to any style imaging table at the front or back. They provide protection against scatter radiation for your lower body and are specially designed to not obstruct positioning of the C-arm.

	22"x24"	22"x26"	22"x28"	22"x32"	22"x36"	24"x24"	24"x36"
Single sided	TS/2224FP	TS/2226FP	TS/2228FP	TS/2232FP	TS/2236FP	TS/2424FP	TS/2436FP